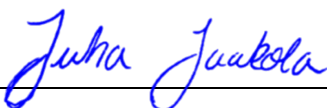



300632

HY-005 Type label specification

Document ID	300632
Version	02

Reviewed by	
Date	2024-05-16
By	Juha Jaakola / Supply Chain Engineer

Approved by	
Date:	2024-05-16
By	Thomas Matzen, QA/RA

The information contained in this document is of a proprietary and confidential nature and is intended only for the persons whom it is directly transmitted by C-RAD. Any reproduction of this document in whole or part, or the divulgence of any of its contents without the prior written consent by C-RAD is prohibited. This document is subject to change without notice.


Document History

Version	Date	Comment
01	2022-02-09	Based on 300284.
02	2024-05-14	HP: Serial number for components and devices of Catalyst+ now indicates if its produced by Innokas or Scanfil. Remove CE and MD mark from Labels

HY-005 Type label

Label specification

The label shall contain the following symbols and text:

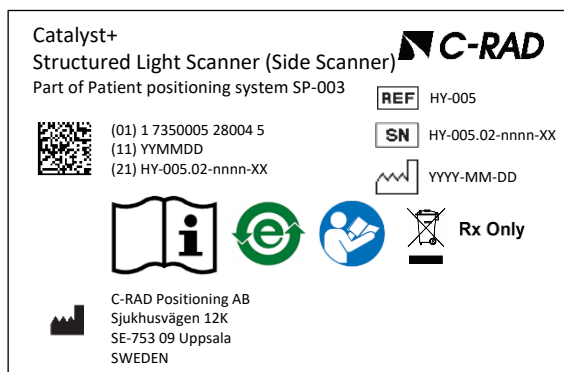

- Update the Manufacturing date “YYMMDD” and “YYYY-MM-DD”, System serial number “nnnn”.
- Update “-XX” with “-SF” only if produced by Scanfil. Otherwise leave blank.
- Material type: B-423 (Brady) white glossy (white polyester).
- Label width: 76,2 mm; height: 50,8 mm
- For every hundred label delivered take a random selection of two labels, attach them to the enclosure and perform the following tests:
 - The label shall be removable with a tool only or by appreciable force.
 - Rub markings by hand, without pressure, first for 15s with a cloth rag soaked with distilled water, then for 15 s with a cloth rag soaked with methylated spirit at ambient temperature and then for 15 s with a cloth rag soaked with isopropyl alcohol. Markings shall be clearly legible after the test.